

## DEVELOPMENT OF A SAMPLE PREPARATION SYSTEM FOR DIOXINS: APPLICATION TO IMMUNOASSAY

Hiroyuki Fujita<sup>1</sup>, Noriaki Hamada<sup>2</sup>, Kazuyuki Sawadaishi<sup>3</sup>, Katsuhisa Honda<sup>1</sup>

<sup>1</sup>Environmental Science for Industry, Ehime University, 3-5-7 Tarumi, Matsuyama Ehime, 790-8566, Japan

<sup>2</sup>Miura Institute of Environmental Science, MIURA Co., Ltd., 864-1, Tsuji, Hojo Ehime, 799-2430, Japan

<sup>3</sup>Kyoto Electronics Manufacturing Co., Ltd., 68 Ninodan-cho, Shinden, Kisshoin, Minami-ku, Kyoto, 601-8317 Japan

### Introduction

Dioxins are highly toxic and bioaccumulative compounds among persistent organic pollutants, and their environmental monitoring and screening are essential. In Japan, conventional testing for dioxins uses the Japan Industrial Standard method (JIS K 0311). This method is expensive and has a relatively slow analytical cycle time. This constrains the number of samples that can be economically tested in projects, leading to uncertainties in the application and interpretation of results for management of site remediation. Immunoassay and bioassay are one way of overcoming these problems<sup>1)</sup>. However, there still remain several problems in clean-up for dioxins contaminated environmental matrices and in analytical complications such as solvent concentration and substitution. For a wide recognition of these assays as monitoring and screening, a suitable pre-treatment method for the assay systems should be developed for all sorts of measurements of analyte.

This study was carried out to develop the automated sample preparation system of dioxins for immunoassay, and the following items were examined; (1) the efficiency of purification by heating multi-layer silica gel column, (2) solvent substitution method using solid-phase adsorbents, and (3) application to immunoassay.

### Methods and Materials

The following study was performed by the experimental instrument, as shown in Fig.1.

#### Samples

Two reference samples of flue gas and fly ash were used for this experiment. Flue gas sample was a mixture of n-Hexane solutions extracted from approx. 30 kinds of flue gas collecting by DioANA<sup>®</sup> filter or an impinger, according to JIS K 0311: 1999. Another fly ash sample was obtained from one type of incinerator. After HCl treatment, dioxins were extracted by soxhlet extraction method and then followed by n-Hexane substitution.

**Purification column (Multi-layer silica gel column) chromatography**


For a validation of the purification column, multi-layer silica gel column was prepared as follows: The packing materials; a few  $\text{Na}_2\text{SO}_4$ (anhydride), 2% KOH silica(0.8g), silica gel(0.05g), 44%  $\text{H}_2\text{SO}_4$  silica(7.3g), silica gel(0.05g), 10%  $\text{AgNO}_3$  silica(3.5g), and a few  $\text{Na}_2\text{SO}_4$ (anhydride), were laminated in sequence into the LC glass column with closed ends (10mm×150mm). Each five ml of one blank and two reference samples was applied on this column, and then eluted with n-Hexane at 60°C or a room temperature.

Each of the eluted solutions was concentrated to 4ml and 100μl using a nitrogen evaporative concentrator and then measured by ultraviolet-visible (V-550 type, Nippon Bunko) and gas chromatography-mass spectrometry (JMS-AUMS200 type, JEOL) respectively.


**Concentration, Separation and Solvent Substitution / Application to Immunoassay**

For investigations of concentration, separation and substitution, an alumina column was set under multi-layer silica gel column (Fig.1). Each reference sample of flue gas and fly ash was applied on a top of multi-layer silica gel column. After heating the column at 60°C, dioxins were eluted with 90ml of n-Hexane kept at 60°C, and then followed by drying the alumina column with  $\text{N}_2$  gas. Thereafter, dioxins in alumina column were reversibly eluted with 1ml of DMSO solution kept at 60°C, and determined by HRGC-HRMS (JMS-700D, JEOL) after n-Hexane substitution. These values of dioxins were compared with the results obtained in conventional treatment, i.e., manual multi-layer silica gel column according to the JIS.

Furthermore, application to immunoassay using flow immunosensor (EndoBioSensor™, Sapidyne Instruments Inc.) and monoclonal anti-dioxin antibody<sup>2,3)</sup> was examined on recovery tests of 2,3,4,7,8-PeCDF added in DMSO solution obtained from the preparation system mentioned above.

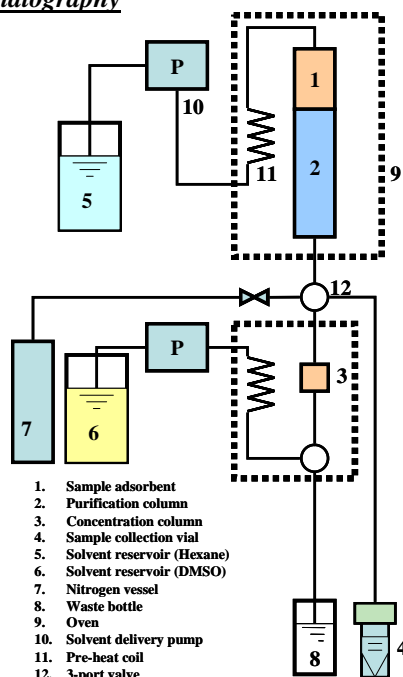

**Fig.1:** Diagram of experimental purification and solvent substitution system.

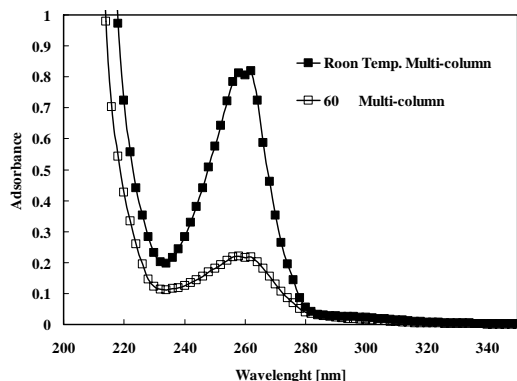

**Fig.2:** UV spectra of flue gas samples purified by multi-layer silica gel column at a room temperature and 60°C.

aromatic and unsaturated compounds and sources of electron transition. While UV spectra at 60°C showed relatively low absorbance, indicating that interferences of immunoassay such as poly aromatic hydrocarbon (PAH) were removed from the flue gas samples. More detailed information was found in the result of GC-MS spectrum (Fig.3). The peaks of PAH, alkyl benzene (ABz) and alkyl biphenyl markedly disappeared. Similar results were also found in fly ash sample.

These observations indicate that progressing in chemical reaction by means of heating multi-layer silica gel column is useful for purification of dioxins analysis.

#### **Concentration, Separation and Solvent Substitution Method**


Recovery tests of dioxins were examined on the automated pre-treatment system (Fig.1) as shown in Fig.4, and their results were compared with those of manual method (multi-layer silica gel column chromatography) according to the JIS (Table 1). High recoveries of dioxins were obtained in this system, as same as those of the JIS method, and their RSD% exhibited smaller values compared to those of the JIS. This indicates that the automated pre-treatment system has a high confidence for dioxins analysis, and also that a process of concentration, separation and solvent substitution using alumina column gives improvement in efficiency.

## **Results and Discussion**


### **Purification column (Multi-layer silica gel column) chromatography**

In multi-layer silica gel column (Fig.1), dioxins were completely eluted with 130ml of n-hexane at a room temperature and with 90ml of n-hexane at 60°C. This indicates that it is possibility to make a short time of analytical cycle and reduce a volume of organic solvent.

UV spectra of flue gas samples purified by multi-layer silicagel column at a room temperature (Fig.2) showed relatively high absorbance of 260nm (230-300nm) and it means presences of

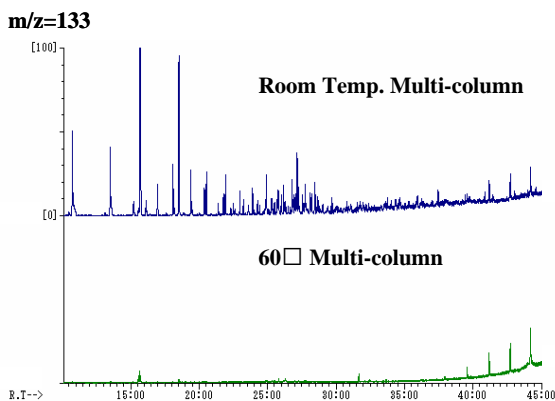


**Fig.3:** GC-LRMS chromatograms of flue gas samples purified by multi-layer silica gel column.

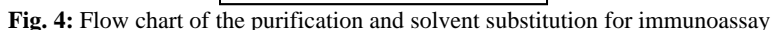

Immunoassay using EndoBioSensor™ was tested on DMSO solutions obtained by the automated system, and their results are shown in Fig.5. The values of 2,3,4,7,8-PeCDF added in DMSO solution linearly increased with an increase in the concentration of 2,3,4,7,8-PeCDF, and recoveries of 2,3,4,7,8-PeCDF were quantitative. The results indicate that the interferences of this immunoassay were removed, further this pre-treatment system is useful for immunoassay.

	Manual (n=3)		Alumina (n=3)	
	Recovery (%)	RSD (%)	Recovery (%)	RSD (%)
2,3,7,8-TeCDD	79	12	93	2
1,2,3,7,8-PeCDD	91	12	100	4
1,2,3,4,7,8-HxCDD	99	14	99	4
1,2,3,6,7,8-HxCDD	104	13	92	3
1,2,3,7,8,9-HxCDD	97	17	106	5
1,2,3,4,6,7,8-HpCDD	122	5	105	4
OCDD	122	2	102	8
2,3,7,8-TeCDF	78	15	92	2
1,2,3,7,8-PeCDF	90	13	103	2
2,3,4,7,8-PeCDF	88	13	90	2
1,2,3,4,7,8-HxCDF	93	15	104	2
1,2,3,6,7,8-HxCDF	100	15	97	3
1,2,3,7,8,9-HxCDF	91	19	87	4
2,3,4,6,7,8-HxCDF	99	13	98	3
1,2,3,4,6,7,8-HpCDF	121	6	102	2
1,2,3,4,7,8,9-HpCDF	113	7	107	3
OCDF	118	1	101	5

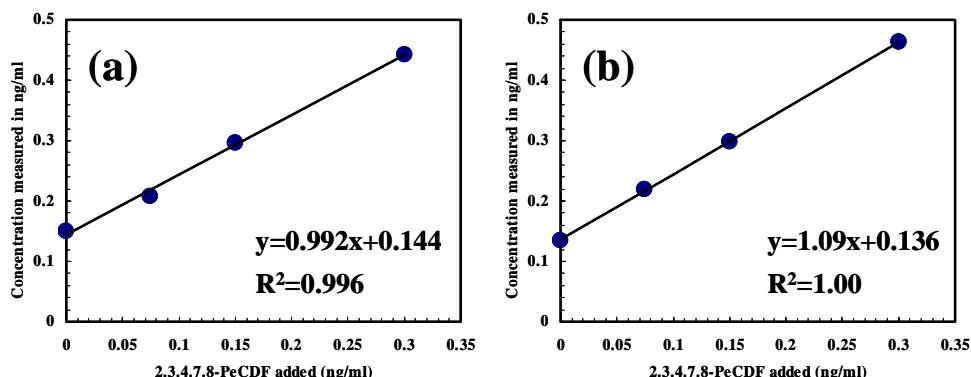

**Fig.5:** Results of recoveries of 2,3,4,7,8-PeCDF from flue gas (a) and fly ash (b) sample. Three different concentrations of 2,3,4,7,8-PeCDF were spiked to 1/10 diluted DMSO solutions prepared on the automated preparation system, and each value of dioxins was measured by EndoBioSensor™.

## Conclusions

An automated pre-treatment system, which is composed of two columns, a multi-layer silica gel column and alumina, is useful for immunoassay of dioxins using EndoBioSensor™. Especially, two special techniques, i.e. heating multi-layer silica gel column and solvent substitution by alumina column, could remove interferences of immunoassay such as aromatic hydrocarbons, and concentrate directly to a small volume of DMSO solution.

## Acknowledgements

This work was partially supported by METI/New Energy and Industrial Technology Development Organization (NEDO).

## References

- 1) Harrison R.O. and Eduljee G. H. (1999) *The Science of the Total Environment* 239, 1-18.
- 2) Glass T.R., Saiki H., Joh T., Taemi Y., Ohmura N. and Lackie S.J. (2004) *Biosensors and Bioelectronics* (in press)
- 3) Takagi Y., Kataoka C., Sawadaishi K., and Shimizu T. (2003) *YAKUGAKU ZASSHI* 123, suppl.1, 84-85 (in Japanese)